

# RapidDoc PinUnlock App

---

November 2022


















# Installing The MTI Pin Unlock App



If necessary, temporarily unplug the scanner to free up a USB-A port

Insert provided USB stick into an available USB-A port

Open the USB Drive Folder and drag and drop the provided ZIP file to a suitable location on the tablet. Unzip the file.

 api-ms-win-crt-string-l1-1-0.dll		12/13/2017 11:54 PM	Application exten...	25 KB
 api-ms-win-crt-time-l1-1-0.dll		12/13/2017 11:54 PM	Application exten...	21 KB
 api-ms-win-crt-utility-l1-1-0.dll		12/13/2017 11:54 PM	Application exten...	19 KB
 base_library		1/12/2018 1:36 PM	Compressed (zipp...	722 KB
 PinUnlock		1/12/2018 1:36 PM	Application	1,582 KB
<input type="checkbox"/> PinUnlock.exe.manifest		1/12/2018 1:36 PM	MANIFEST File	2 KB
 pyexpat.pyd		12/13/2017 11:54 PM	PYD File	186 KB
 python36.dll		12/13/2017 11:54 PM	Application exten...	3,519 KB


The MTI Pin Unlock App is this icon

Pin the Pin Unlock icon to the Start Menu and/or the Task Bar



# MTI Pin Unlock App

- Default Pin Code is 97124
- Default Admin Pin is 0000
- To undock the tablet, type 97124, then press Unlock
- If a proper PIN is entered, a green indicator will appear in the App and the RapidDoc VESA mount will open the latches with a clear clicking sound

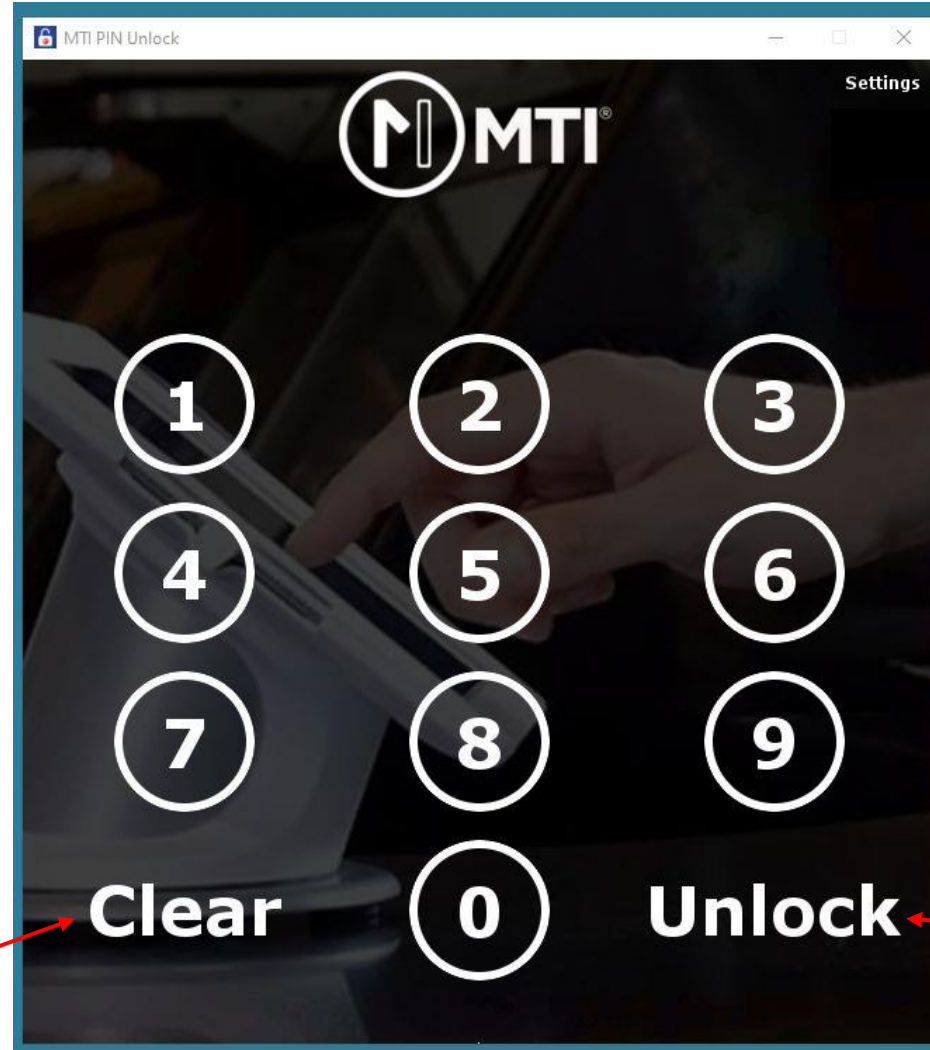


Settings  
button used  
to enter  
Admin Mode


Type a PIN,  
then press  
Unlock

Press Clear  
to start over

# MTI Pin Unlock App – Adding Users

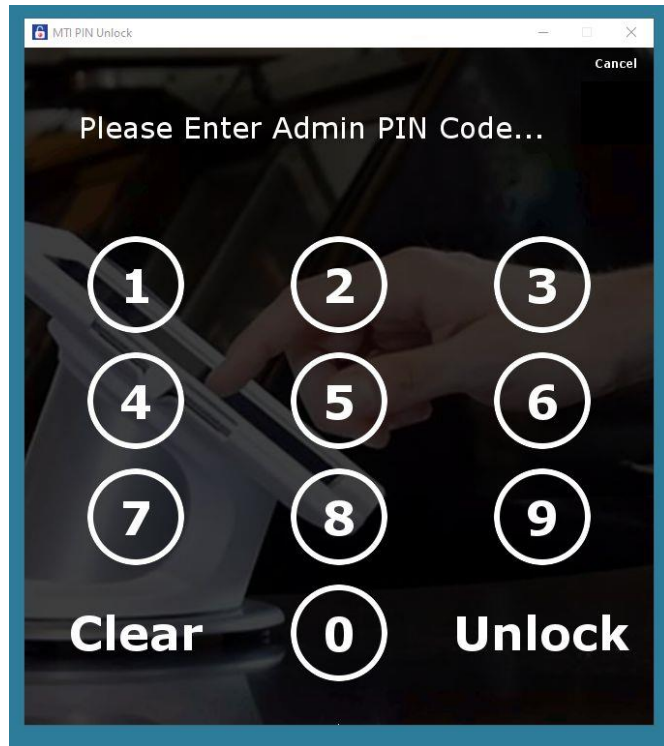


Click on the “Settings” button in the top right corner to enter Admin Mode  
Type “0000” then press Unlock

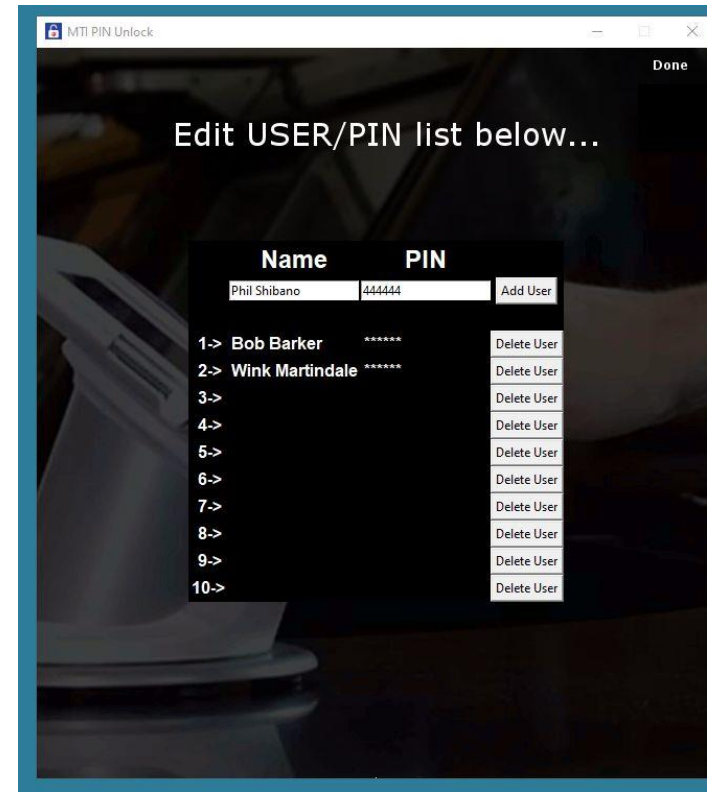

In Admin Mode, enter Names and New Pins as desired. Delete old users in this screen too. Hit “Done” button in top right corner when complete.  
PIN numbers are stored locally on tablet.



# Undocking the Tablet with the PinUnlock App

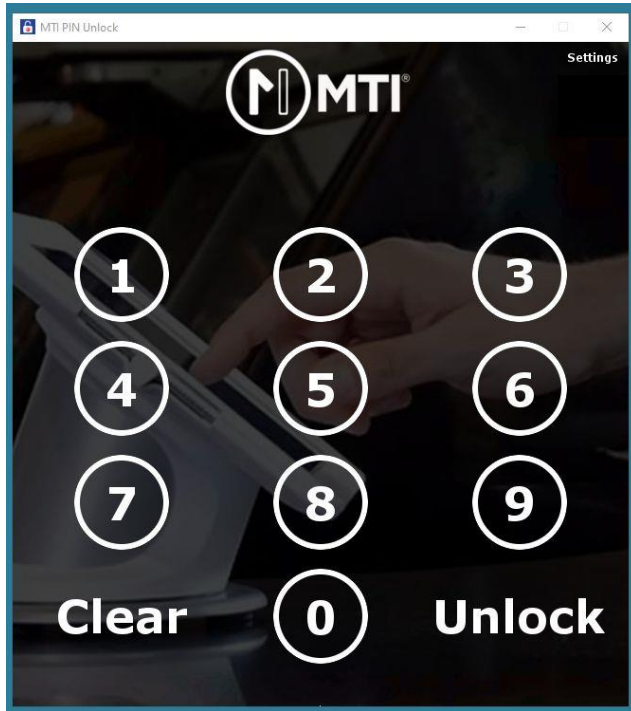

Type a stored PIN number or the default PIN and push the Unlock button.

If the PIN is accepted, Green dots will be shown  
If the PIN is not accepted Red dots will be shown  
You should hear a “click” from the RapidDoc Mount as the latches disengage. The latches will stay open for 5 seconds, then re-close automatically



Peel Right to Left  
To Release Tablet

Even though the tablet is unlatched, magnets prevent it from falling. To remove the tablet, simply “peel” the tablet from left to right or right to left to break the magnetic force



To redock, simply align the tablet mount with the RapidDoc Mount. It will magnetically guide into place and the latches will capture the tablet automatically



*AMERICA'S PRODUCTIVITY PARTNER*

*gbscorp.com | EasyID Support: 800-369-0595*

